

# Mercodia

# Ovine Insulin

# ELISA

Directions for Use

**10-1202-01**

**Reagents for 96 determinations**





For Research Use Only

Manufactured by

**Mercodia AB**

Sylveniusgatan 8A  
SE-754 50 Uppsala  
Sweden

## Explanation of symbols used on labels

 $\Sigma = 96$	Reagents for 96 determinations
	Expiry date
	Store between 2–8°C
	Lot No.

## **Intended Use**

Mercodia Ovine Insulin ELISA provides a method for the quantitative determination of insulin in ovine serum and plasma.

## **Summary and explanation of the test**

Insulin is the principal hormone responsible for the control of glucose metabolism. It is synthesised in the  $\beta$ -cells of the islets of Langerhans as the precursor, proinsulin, which is processed to form C-peptide and insulin. Both are secreted in equimolar amounts into the portal circulation. The mature insulin molecule comprises two polypeptide chains, the A chain and the B chain. The two chains are linked together by two inter-chain disulphide bridges. There is also an intra-chain disulphide bridge in the A chain.

Secretion of insulin is mainly controlled by plasma glucose concentration, and the hormone has a number of important metabolic actions. Its principal function is to control the uptake and utilisation of glucose in peripheral tissues via the glucose transporter. This and other hypoglycaemic activities, such as the inhibition of hepatic gluconeogenesis and glycogenolysis are counteracted by the hyperglycaemic hormones including glucagon, epinephrine (adrenaline), growth hormone and cortisol.

Extensive research on how to improve the nutritional, metabolic and health status of ruminants has been at focus for a long time. A decrease in dry matter intake is a major physiological change in ruminants. This may lead to several metabolic disorders such as ketosis, fatty liver and hypocalcemia<sup>1,2</sup>. Food intake is a complex mechanism, regulated by several factors including for example hormones, metabolites and environmental factors<sup>1-3</sup>.

## **Principle of the procedure**

Mercodia Ovine Insulin ELISA is a solid phase two-site enzyme immunoassay. It is based on the direct sandwich technique in which two monoclonal antibodies are directed against separate antigenic determinants on the insulin molecule. During incubation, insulin in the sample reacts with peroxidase-conjugated anti-insulin antibodies and anti-insulin antibodies bound to the microtitration well. After a simple washing step that removes unbound enzyme labelled antibody, the bound conjugate is detected by reaction with 3,3',5,5'-tetramethylbenzidine (TMB). The reaction is stopped by the addition of acid, giving a colorimetric endpoint that can be read spectrophotometrically.

## Warnings and precautions

- All samples should be handled as capable of transmitting infections.
- Not for internal or external use in humans or animals.
- Each well can only be used once.
- For Research Use Only.
- The Stop Solution contains <5% Sulphuric acid.

The Stop Solution is labeled:



### Danger

H318 – Causes serious eye damage.

H315 – Causes skin irritation.

P280 – Wear protective gloves. Wear eye or face protection.

P264 – Wash hands thoroughly after handling.

P302 + P352 + P362 + P364 – IF ON SKIN: Wash with plenty of soap and water. Take off contaminated clothing and wash it before reuse.

P332 + P313 – If skin irritation occurs: Get medical attention.

P305 + P351 + P338 + P310 – IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a POISON CENTER or physician.

- The Enzyme Conjugate Buffer, Cal 0, 1, 2, 3, 4, 5 and Wash Buffer contain <0.06% 5-chloro-2-methyl-4-isothiazolin-3-one and 2-methyl-2H-isothiazol-3-one (3:1).

The Enzyme Conjugate Buffer, the Calibrators and Wash Buffer are labeled:



### Warning

H317 – May cause an allergic skin reaction.

P280 – Wear protective gloves.

P261 – Avoid breathing vapour.

P272 – Contaminated work clothing should not be allowed out of the workplace.

P302 + P352 – IF ON SKIN: Wash with plenty of soap and water.

P333 + P313 – If skin irritation or rash occurs: Get medical attention.

P501 – Dispose of contents and container in accordance with all local, regional, national and international regulations.

## Material required but not provided

- Pipettes with appropriate volumes (repeating pipettes preferred for addition of enzyme conjugate 1X solution, Substrate TMB and Stop Solution)
- Tubes, beakers and cylinders for reagent preparation
- Redistilled water
- Magnetic stirrer
- Vortex mixer
- Microplate reader with 450 nm filter
- Microplate shaker (700–900 cycles per minute, orbital movement)
- Microplate washing device with overflow function (recommended but not required)

## Reagents

Each Mercodia Ovine Insulin ELISA kit contains reagents for 96 wells, sufficient for 42 samples and one calibrator curve in duplicate. For larger series of assays, use pooled reagents from packages bearing identical lot numbers. The expiry date for the complete kit is stated on the outer label. The recommended storage temperature is 2–8°C.

<b>Coated Plate</b> Mouse monoclonal anti-insulin For unused microtitration strips, reseal the bag using adhesive tape, store at 2–8°C and use within 2 months.	1 plate	96 wells 8-well strips	Ready for use
<b>Calibrators 1, 2, 3, 4, 5</b> Insulin Concentration stated on vial label. Color coded yellow	5 vials	1000 µL	Ready for Use
<b>Calibrator 0</b> Color coded yellow	1 vial	5 mL	Ready for use
<b>Enzyme Conjugate 11X</b> Peroxidase conjugated mouse monoclonal anti-insulin	1 vial	1.3 mL	Preparation, see below
<b>Enzyme Conjugate Buffer</b> Color coded blue	1 vial	13 mL	Ready for use
<b>Wash Buffer 21X</b> Storage after dilution: 2–8°C for 2 months	1 bottle	50 mL	Dilute with 1000 mL redistilled water to make wash buffer 1X solution.
<b>Substrate TMB</b> Colorless solution <i>Note! Light sensitive!</i>	1 bottle	22 mL	Ready for use
<b>Stop Solution</b> 0.5 M H <sub>2</sub> SO <sub>4</sub>	1 vial	7 mL	Ready for use

## Preparation of enzyme conjugate 1X solution

Prepare the needed volume of enzyme conjugate 1X solution by dilution of Enzyme Conjugate 11X, (1+10) in Enzyme Conjugate Buffer according to the table below. Mix gently. When preparing enzyme conjugate 1X solution for the whole plate or if the reagents are to be used within 2 weeks, pour all of the Enzyme Conjugate Buffer into the Enzyme Conjugate 11X vial.

Number of strips	Enzyme Conjugate 11X	Enzyme Conjugate Buffer
12 strips	1 vial	1 vial
8 strips	700 µL	7 mL
4 strips	400 µL	4 mL

Storage after dilution: 2–8°C. Use within 2 weeks.

## Specimen collection and handling

### Serum

Collect blood by venipuncture, allow to clot and separate the serum by centrifugation. Samples can be stored at 2–8°C up to 24 hours. For longer periods, store samples at -20°C. Avoid repeated freezing and thawing.

### Plasma

Collect blood by venipuncture into tubes containing heparin, citrate or EDTA as anticoagulant, and separate the plasma fraction by centrifugation. Samples can be stored at 2–8°C up to 24 hours. For longer periods store samples at -20°C. Avoid repeated freezing and thawing.

## Preparation of samples

No dilution is normally required for serum or plasma. All samples containing ovine insulin above the highest calibrator should be diluted with Calibrator O or with Merckodia Diabetes Sample Buffer 10-1195-01.

## Test procedure

All reagents and samples must be brought to room temperature before use. Perform each determination in duplicate for calibrators and samples. Prepare a calibrator curve for each assay run. The product has been optimized and validated without plate sealer.

1. Prepare enzyme conjugate 1X solution (according to the table on previous page) and wash buffer 1X solution.
2. Prepare sufficient microplate wells to accommodate Calibrators, controls and samples in duplicate.
3. Pipette 25  $\mu$ L each of Calibrators, controls and samples into appropriate wells.
4. Add 100  $\mu$ L of enzyme conjugate 1X solution into each well.
5. Incubate on a plate shaker (700–900 rpm) for 2 hours at room temperature (18–25°C).
6. Wash 6 times with 700  $\mu$ L wash buffer 1X solution per well using an automatic plate washer with overflow-wash function, after final wash, invert and tap the plate firmly against absorbent paper. Do not include soak step in washing procedure.  
Or manually,  
Discard the reaction volume by inverting the microplate over a sink.  
Add 350  $\mu$ L wash buffer 1X solution to each well. Discard the wash buffer 1X solution, tap firmly several times against absorbent paper to remove excess liquid. Repeat 5 times. Avoid prolonged soaking during washing procedure.
7. Add 200  $\mu$ L Substrate TMB into each well.
8. Incubate on the bench for 15 minutes at room temperature (18–25°C).
9. Add 50  $\mu$ L Stop Solution to each well.  
Place the plate on the shaker for approximately 5 seconds to ensure mixing.
10. Read optical density at 450 nm and calculate results.  
Read within 30 minutes.

*Note!* Be extra careful not to contaminate the Substrate TMB with enzyme conjugate solution!

## Internal quality control

Commercial control such as Mercodia Insulin Control Animal Low, Medium and High (10-1221-01) and/or internal serum pools with low, intermediate and high ovine insulin concentrations should routinely be assayed as samples, and results charted from day to day. It is good laboratory practice to record the following data for each assay: kit lot number; dilution and/or reconstitution dates of kit components; OD values for the blank, Calibrators and controls.

## Calculation of results

The concentration of ovine insulin is obtained by computerized data reduction of the absorbance for the Calibrators, except Calibrator 0, versus the concentration using cubic spline regression.


## Example of results

Wells	Identity	A <sub>450 nm</sub>	Mean conc. µg/L
1A-B	Calibrator 0	0.058/0.059	
1C-D	Calibrator 1*	0.093/0.090	
1E-F	Calibrator 2*	0.157/0.155	
1G-H	Calibrator 3*	0.412/0.414	
2A-B	Calibrator 4*	1.301/1.309	
2C-D	Calibrator 5*	2.765/2.735	
2E-F	Sample 1	0.193/0.190	0.206
2G-H	Sample 2	0.356/0.361	0.437
3A-B	Sample 3	0.788/0.764	0.921

\*Exact concentration stated on vial label.

## Example of calibrator curve

A typical calibrator curve is shown here. Do not use this curve to determine actual assay results.

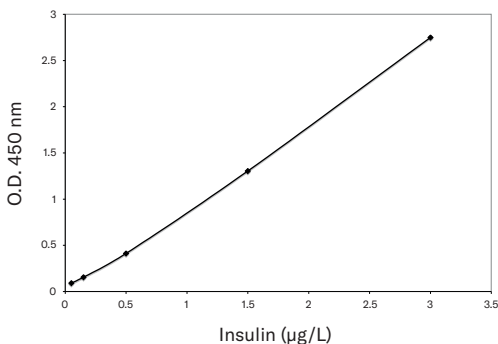



## Limitations of the procedure

Grossly lipemic, icteric or hemolyzed samples do not interfere in the assay. Insulin is, however, degraded over time in hemolyzed samples. The degradation could give falsely low values and contributes to higher inter assay variation.

## Expected values

Good practice dictates that each laboratory establishes its own expected range of values.

## Performance characteristics

### Detection limit

Detection limit is defined as the Capability of Detection according to ISO11843-Part 1. Capability of Detection should be seen as part of a method validation, rather than the lowest concentration that can be measured.

The detection limit is 0.025 ( $\mu\text{g/L}$ ) determined with the methodology described in ISO11843-Part 4.

Concentration of samples with absorbance below Calibrator 1 should not be calculated, instead expressed as less or equal to  $\leq$  the concentration indicated on the vial for Calibrator 1.

### Recovery

Recovery upon addition is 94% – 114% (mean 103%).

Recovery upon dilution is 68% – 122% (mean 89%).

### Hook effect

Samples with a concentration of up to 1 000  $\mu\text{g/L}$  can be measured without giving falsely low results.

### Precision

Each sample was analyzed in 4 replicates on 19 different occasions.

Sample	Mean value $\mu\text{g/L}$	Coefficient of variation	
		Repeatability %*	Within laboratory %**
1	0.201	3.7	6.8
2	0.403	1.2	4.6
3	0.859	1.7	4.8

\*Within assay variation

\*\*Total assay variation

## Specificity

The following cross reactions have been found:

	Crossreaction
Human insulin	133%
NovoRapid	2.0%
Lantus	21%
Humalog	n.d.
Levemir	n.d.

n.d. = not detected

## Calibration

Mercodia Ovine Insulin ELISA is calibrated against an in house reference preparation of ovine insulin.

## Warranty

The performance data presented here was obtained using the procedure indicated. Any change or modification in the procedure not recommended by Mercodia AB may affect the results, in which event Mercodia AB disclaims all warranties expressed, implied or statutory, including the implied warranty of merchantability and fitness for use.

Mercodia AB and its authorized distributors, in such event, shall not be liable for damages indirect of consequential.

## References

1. Hagino A, Inomata E, Sato T, Ohtomo Y, Sasaki Y, Obara Y (2005) Effect in sheep of dietary concentrate content on secretion of growth hormone, insulin and insulin-like growth factor-1 after feeding. *Ani Sci J* 76:55-63
2. McCann JP, Loo SC, Aalseth DL, Aribat T (1997) Differential effects of GH stimulation on fasting and prandial metabolism and plasma IGFs and IGF-binding proteins in lean and obese sheep. *J Endocrin* 154:329-346
3. Melandez p, Krueger T, White J, Badinga L, Verstegen J, donovan GA, Archbald LF (2006) Effect of ghrelin in dry matter intake and energy metabolism in parturient sheep: a preliminary study. *Theriogenology* 66: 1961-1968

Further references can be found on our website: [www.mercodia.com](http://www.mercodia.com)

## Notes

[illegible]

## Notes

This image shows a single page of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

[illegible]

## Notes

[illegible]

	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C												
D												
E												
F												
G												
H												

**Experiment:**

**Date:**

**KIT LOT#:**

**Summary of protocol sheet**  
**Mercodia Ovine Insulin ELISA**

Add Calibrators, controls* and Samples	25 µL
Add enzyme conjugate 1X solution	100 µL
Incubate	2 hours at 18–25°C on a plate shaker (700–900 rpm)
Wash plate with wash buffer 1X solution	700µL, 6 times
Add Substrate TMB	200 µL
Incubate	15 minutes at 18–25°C
Add Stop Solution	50 µL Shake for 5 seconds to ensure mixing
Measure A <sub>450 nm</sub>	Evaluate results

\*not included

For full details see page 7

For technical support please contact: [support@merckodia.com](mailto:support@merckodia.com)